

**BioVendor**

Research  
and Diagnostic Products



## Human CTRP5 ELISA

Product Data Sheet

Cat. No.: RAG003R

For Research Use Only

## CONTENTS

---

1.	INTENDED USE	3
2.	HANDLING, STORAGE	3
3.	INTRODUCTION	3
4.	TEST PRINCIPLE	4
5.	TECHNICAL HINTS	4
6.	REAGENT SUPPLIED	5
7.	MATERIALS REQUIRED BUT NOT SUPPLIED	5
8.	PREPARATION OF REAGENTS	6
9.	PREPARATION OF SAMPLES	7
10.	ASSAY PROCEDURE	8
11.	CALCULATIONS	9
12.	PERFORMANCE CHARACTERISTICS	10
13.	TROUBLESHOOTING	12
14.	REFERENCES	13

➤➤ **This kit is manufactured by:**  
**BioVendor – Laboratorní medicína a.s.**

➤➤ **Use only the current version of Product Data Sheet enclosed with the kit!**

## 1. INTENDED USE

---

The CTRP5 (human) Competitive ELISA Kit is to be used for the *in vitro* quantitative determination of human CTRP5 in serum, plasma and cell culture supernatant. This ELISA Kit is for research use only.

## 2. HANDLING, STORAGE

---

- Reagent must be stored at 2-8°C when not in use.
- Plate and reagents should be at room temperature before use.
- Do not expose reagents to temperatures greater than 25°C.

## 3. INTRODUCTION

---

CTRP5 is a member of the C1q and tumor necrosis factor superfamily whose structure resembles adiponectin, being a 25 kD secretory protein (1-2). CTRP5 was identified in a cDNA library enriched for genes showing expression specific for the retinal pigment epithelium (RPE) (3). A subsequent RT-PCR study demonstrated CTRP5 expression in RPE, liver, lung, placenta, and brain. It has been proposed that a CTRP mutation (S163R) plays a critical role in affecting its higher order protein structure, potentially leading to a cause of abnormal adhesion between the RPE and Bruch membrane (4-5). A recent data showed that CTRP is one of the genes highly induced by elimination of mitochondria and was able to activate AMPK in a rat myotube cell line, L6. Stimulation of L6 with recombinant CTRP5, full length as well as globular domain, enhanced glucose uptake and fatty acid oxidation. These biochemical events did not seem to be mediated via AdipoR1 or AdipoR2, suggesting that a novel receptor(s) may exist for CTRP5 in this muscle cell line (6). Some CTRP members can physically interact with adiponectin, forming various multimeric structures (2). A compensatory upregulation of mouse CTRP members was observed in adiponectin-deficient mice (1). Therefore, measuring serum or plasma CTRP5 may provide important information on its involvement in novel metabolism.

## 4. TEST PRINCIPLE

---

This assay is a competitive Enzyme Linked-Immunosorbent Assay (ELISA) for quantitative determination of human CTRP5 in biological fluids. A polyclonal antibody recognizing native human CTRP5 reacts with a series of predetermined recombinant human CTRP5 standard proteins or samples under competition in the human CTRP5-coated plate. Their relative reactivity is plotted with that of the standard proteins.

## 5. TECHNICAL HINTS

---

- It is recommended that all standards, QC samples and samples be run in duplicate.
- Do not combine leftover reagents with those reserved for additional wells.
- Reagents from the kit with a volume less than 100 µl should be centrifuged.
- Residual wash liquid should be drained from the wells after last wash by tapping the plate on absorbent paper.
- Crystals could appear in the 10X solution due to high salt concentration in the stock solutions. Crystals are readily dissolved at room temperature or at 37°C before dilution of the buffer solutions.
- Once reagents have been added to the 8-well strips, DO NOT let the strips DRY at any time during the assay.
- Keep TMB Substrate Solution protected from light.
- The Stop Solution consists of phosphoric acid. Although diluted, the Stop Solution should be handled with gloves, eye protection and protective clothing.

## 6. REAGENT SUPPLIED

<b>Kit Components</b>	<b>Quantity</b>
1 plate coated with human CTRP5 Recombinant Protein	12 x 8-well strips
1 bottle Wash Buffer 10X	50 ml
1 bottle Diluent 5X	50 ml
1 bottle Detection Antibody	12 ml
1 vial Detector 100X (HRP Conjugated anti-rabbit IgG)	150 $\mu$ l
1 vial human CTRP5 Standard (lyophilized)	5 $\mu$ g
1 vial human CTRP5 QC sample (lyophilized)	
1 bottle TMB Substrate Solution	12 ml
1 bottle Stop Solution	12 ml
3 plate sealers (plastic film)	

## 7. MATERIALS REQUIRED BUT NOT SUPPLIED

- Microtiterplate reader at 450 nm, with the correction wavelength set at 540 nm or 570 nm
- Calibrated precision single and multi-channel pipettes. Disposable pipette tips
- Deionized water
- Microtubes or equivalent for preparing dilutions
- Disposable plastic containers for preparing working buffers
- Plate washer: automated or manual
- Glass or plastic tubes for diluting and aliquoting standard

## 8. PREPARATION OF REAGENTS

---

**NOTE:** Prepare just the appropriate amount of the buffers necessary for the assay.

- **Wash Buffer 10X** has to be diluted with deionized water 1:10 before use (e.g. 50 ml Wash Buffer 10X + 450 ml water) to obtain Wash Buffer 1X.
- **Diluent 5X** has to be diluted with deionized water 1:5 before use (e.g. 50 ml Diluent 5X + 200 ml water) to obtain Diluent 1X.
- **Detector 100X (HRP Conjugated anti-rabbit IgG)** has to be diluted to the working concentration by adding 120 µl in 12 ml of Diluent 1X (1:100).

**NOTE:** The diluted Detector is used within one hour of preparation.


- **Human CTRP5 Standard (STD)** has to be reconstituted with 1 ml of deionized water.
  - This reconstitution produces a stock solution of 5 µg/ml. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes. Mix well prior to making dilutions.

**NOTE:** The reconstituted standard is aliquoted and stored at -20°C.

- Dilute the standard protein concentrate (STD) (**5 µg/ml**) in Diluent 1X. A seven-point standard curve in Diluent 1X is recommended.
  - Suggested standard points are:  
**5 , 2.5 , 1 , 0.5 , 0.25 , 0.1 , 0.01 and 0.001 µg/ml.**
- **Human CTRP5 QC sample** has to be reconstituted with 1 ml of deionized water.
    - Refer to the Certificate of Analysis for current QC sample concentration. Mix the QC sample to ensure complete reconstitution and allow the QC sample to sit for a minimum of 15 minutes. The reconstituted QC sample is ready to use, do not dilute it.

**Dilute further for the standard curve:**

To obtain	Add	Into
5 µg/ml	-	-
2.5 µg/ml	300 µl of CTRP5 (5 µg/ml)	300 µl of Diluent 1X
1 µg/ml	200 µl of CTRP5 (2.5 µg/ml)	300 µl of Diluent 1X
0.5 µg/ml	300 µl of CTRP5 (1 µg/ml)	300 µl of Diluent 1X
0.25 µg/ml	300 µl of CTRP5 (0.5 µg/ml)	300 µl of Diluent 1X
0.1 µg/ml	200 µl of CTRP5 (0.25 µg/ml)	300 µl of Diluent 1X
0.01 µg/ml	50 µl of CTRP5 (0.1 µg/ml)	450 µl of Diluent 1X
0.001 µg/ml	50 µl of CTRP5 (0.01 µg/ml)	450 µl of Diluent 1X

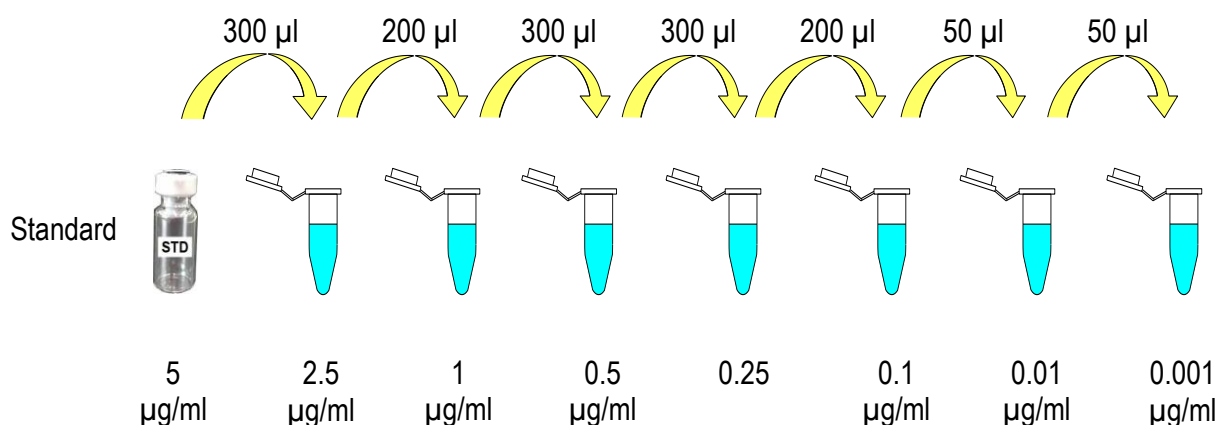

## 9. PREPARATION OF SAMPLES

**Serum** : Use a serum separator tube. Let samples clot at room temperature for 30 minutes before centrifugation for 20 minutes at 1,000xg. Assay freshly prepared serum or store serum in aliquot at ≤ -20°C for later use. Avoid repeated freeze/thaw cycles.

**Plasma** : Collect plasma using heparin, EDTA, or citrate as an anticoagulant. Centrifuge for 15 minutes at 1000xg within 30 minutes of collection. Assay freshly prepared plasma or store plasma sample in aliquot at ≤ -20°C for later use. Avoid repeated freeze/ thaw cycles.

**Serum, Plasma or Cell Culture Supernatant** have to be diluted in Diluent 1X. Samples containing visible precipitates must be clarified before use.

**NOTE:** As a starting point, 1/2 dilution of serum or plasma is recommended! If samples fall the outside range of assay, a lower or higher dilution may be required!

## 10. ASSAY PROCEDURE

---

1. Determine the number of 8-well strips needed for the assay and insert them in the frame for current use. The extra strips should be resealed in the foil pouch bag and stored at 4°C.  
**NOTE:** Remaining 8-well strips coated with CTRP5 protein when opened can be stored at 4°C for up to 1 month.
2. Add 50 µl of the different standards and reconstituted QC sample into the appropriate wells in duplicate! At the same time, add 50 µl of diluted serum, plasma or cell culture supernatant samples in duplicate to the wells (**see 8 Preparation of Reagents and 9 Preparation of Samples**).
3. Add 50 µl to each well of the Detection Antibody and tap gently on the side of the plate to mix
4. Cover the plate with plate sealer and incubate for **1 hour at 37°C**.
5. Aspirate the coated wells and add 300 µl of Wash Buffer 1X using a multi-channel pipette or auto-washer. Repeat the process for a total of three washes. After the last wash, complete removal of liquid is essential for good performance.
6. Add 100 µl to each well of the diluted Detector (**see 8 Preparation of Reagents**).
7. Cover the plate with plate sealer and incubate for **1 hour at 37°C**.
8. Aspirate the coated wells and add 300 µl of Wash Buffer 1X using a multi-channel pipette or auto-washer. Repeat the process for a total of five washes. After the last wash, complete removal of liquid is essential for good performance.
9. Add 100 µl to each well of TMB Substrate Solution.
10. Allow the color reaction to develop **at room temperature (RT°C) in the dark for 10 minutes**
11. Stop the reaction by adding 100 µl of Stop Solution. Tap the plate gently to ensure thorough mixing. The substrate reaction yields a blue solution that turns yellow when Stop Solution is added.

**! CAUTION: CORROSIVE SOLUTION!**


12. Measure the OD at 450 nm in an ELISA reader within 30 minutes.



## 11. CALCULATIONS

- Average the duplicate readings for each standard, QC and sample.
- Generate the standard curve by plotting the average absorbance obtained for each standard concentration on the vertical (Y) axis vs. the corresponding CTRP5 concentration ( $\mu\text{g/ml}$ ) on the horizontal (X) axis.
- Calculate the CTRP5 concentrations of samples by interpolation of the regression curve formula as shown above in a form of a 4-parameter logistic equation.
- If the test samples were diluted, multiply the interpolated values by the dilution factor to calculate the concentration of human CTRP5 in the samples.

The following data are obtained using the different concentrations of standard as described in this protocol:

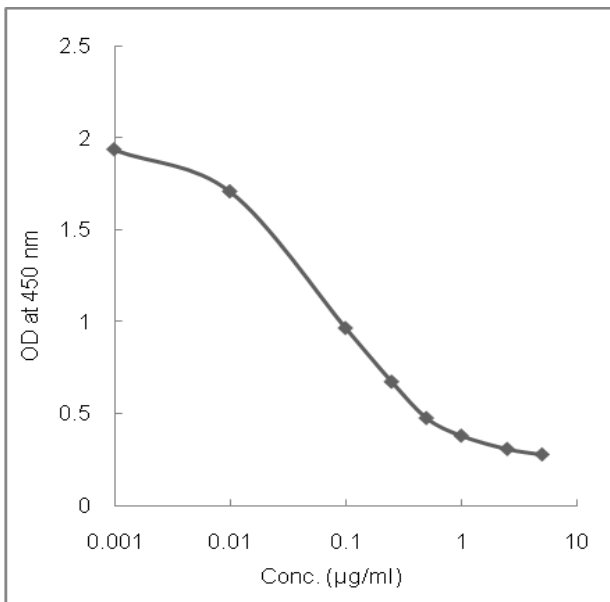

**Figure:** Standard curve

Standard hCTR5 ( $\mu\text{g/ml}$ )	Optical Density (mean)
5	0.277
2.5	0.307
1	0.379
0.5	0.475
0.25	0.672
0.1	0.965
0.01	1.707
0.001	1.935

## 12. PERFORMANCE CHARACTERISTICS

Typical analytical data of BioVendor human CTRP5 ELISA, Clinical Range are presented in this chapter

- **Sensitivity** (Limit of detection)

The lowest level of CTRP5 that can be detected by this assay is 1 ng/ml.

**NOTE:** The Limit of detection was measured by adding two standard deviations to the mean value of 50 zero standard.

- **Assay range**

0.001 µg/ml – 5 µg/ml

- **Specificity**

This ELISA is specific for the measurement of natural and recombinant human CTRP5. It does not cross-react with mouse CTRP5, human CTRP2 (globular), human CTRP6, human CTRP7 (globular), human CTRP9 (globular), human CTRP10 (globular), mouse CTRP2 (globular), mouse CTRP9 (globular), human adiponectin, human adiponectin (globular), mouse adiponectin, mouse adiponectin (globular), rat adiponectin, rat adiponectin (globular), human RBP4, human Nampt, human vaspin, human GPX3, human ANGPTL3, human progranulin, human leptin.

Human CTRP5 (globular) and CTRP6 (globular) show weakly 20% cross-reactivity in this assay.

- **Precision:**

*Intra-assay* (n =9)

Four samples of known concentrations of human CTRP5 were assayed in replicates 9 times to test precision within an assay.

Sample	Mean (ng/ml)	SD (ng/ml)	CV (%)
1	0.351	0.02	5.87
2	0.293	0.03	10.00
3	0.250	0.01	4.05
4	0.165	0.01	6.06

*Inter-assay (n = 3)*

Four samples of known concentrations of human CTRP5 were assayed in 3 separate assays to test precision between assays.

Sample	Mean (ng/ml)	SD (ng/ml)	CV (%)
1	0.348	0.02	6.36
2	0.279	0.03	9.08
3	0.253	0.02	8.11
4	0.162	0.01	6.77

- **Linearity**

Different human serum samples containing CTRP5 were diluted several fold (1 to 1/4) and the measured recoveries ranged from 79% to 130%.

Sample	Dilution	Observed (ng/ml)	Expected (ng/ml)	% of Expected
1	1	0.260	0.260	100
	1 : 2	0.125	0.130	95
	1 : 4	0.071	0.065	108
2	1	0.174	0.174	100
	1 : 2	0.082	0.087	94
	1 : 4	0.056	0.044	127
3	1	0.321	0.321	100
	1 : 2	0.127	0.161	79
	1 : 4	0.071	0.080	88

- **Comparison of serum samples with plasma samples**

Different human serum samples containing CTRP5 were compared with plasma samples.

Sample	Serum (µg/ml)	Plasma (µg/ml)		
		Citrate	EDTA	Heparin
1	0.127	0.374	0.343	0.126
2	0.098	0.196	0.346	0.098
3	0.125	0.21	0.504	0.134
4	0.145	0.31	0.41	0.152

- **Expected values:**

CTRP5 levels range in plasma and serum from **0.05 to > 0.5 µg/ml** (from healthy donors).

### 13. TROUBLESHOOTING

<b>PROBLEM</b>	<b>POSSIBLE CAUSES</b>	<b>SOLUTIONS</b>
No signal or weak signal	Omission of key reagent	Check that all reagents have been added in the correct order.
	Washes too stringent	Use an automated plate washer if possible.
	Incubation times inadequate	Incubation times should be followed as indicated in the manual.
	Plate reader settings not optimal	Verify the wavelength and filter setting in the plate reader.
	Incorrect assay temperature	Use recommended incubation temperature. Bring substrates to room temperature before use.
High background	Concentration of detector too high	Use recommended dilution factor.
	Inadequate washing	Ensure all wells are filling wash buffer and are aspirated completely.
Poor standard curve	Wells not completely aspirated	Completely aspirate wells between steps.
	Reagents poorly mixed	Be sure that reagents are thoroughly mixed.
Unexpected results	Omission of reagents	Be sure that reagents were prepared correctly and added in the correct order.
	Dilution error	Check pipetting technique and double-check calculations.


## 14. REFERENCES

---


### References to CTRP5:

1. Molecular, biochemical and functional characterizations of C1q/TNF family members: adipose-tissue-selective expression patterns, regulation by PPAR-gamma agonist, cysteine-mediated oligomerizations, combinatorial associations and metabolic functions: G.W. Wong, et al.; *Biochem. J.* 416, 161 (2008)
2. Identification and characterization of CTRP9, a novel secreted glycoprotein, from adipose tissue that reduces serum glucose in mice and forms heterotrimers with adiponectin: G.W. Wong, et al.; *FASEB J.* 23, 241 (2009)
3. CTRP5 is a membrane-associated and secretory protein in the RPE and ciliary body and the S163R mutation of CTRP5 impairs its secretion: M.N. Mandal, et al.; *Invest. Ophthalmol. Vis. Sci.* 47, 5505 (2006)
4. Spatial and temporal expression of MFRP and its interaction with CTRP5: M.N. Mandal, et al.; *Invest. Ophthalmol. Vis. Sci.* 47, 5514 (2006)
5. Mutation in a short-chain collagen gene, CTRP5, results in extracellular deposit formation in late-onset retinal degeneration: a genetic model for age-related macular degeneration: C. Hayward, et al.; *Hum. Mol. Genet.* 12, 2657 (2003)
6. C1q tumor necrosis factor alpha-related protein isoform 5 is increased in mitochondrial DNA-depleted myocytes and activates AMP-activated protein kinase: S.Y. Park, et al.; *J. Biol. Chem.* 284, 27780 (2009)

## Assay Procedure Summary

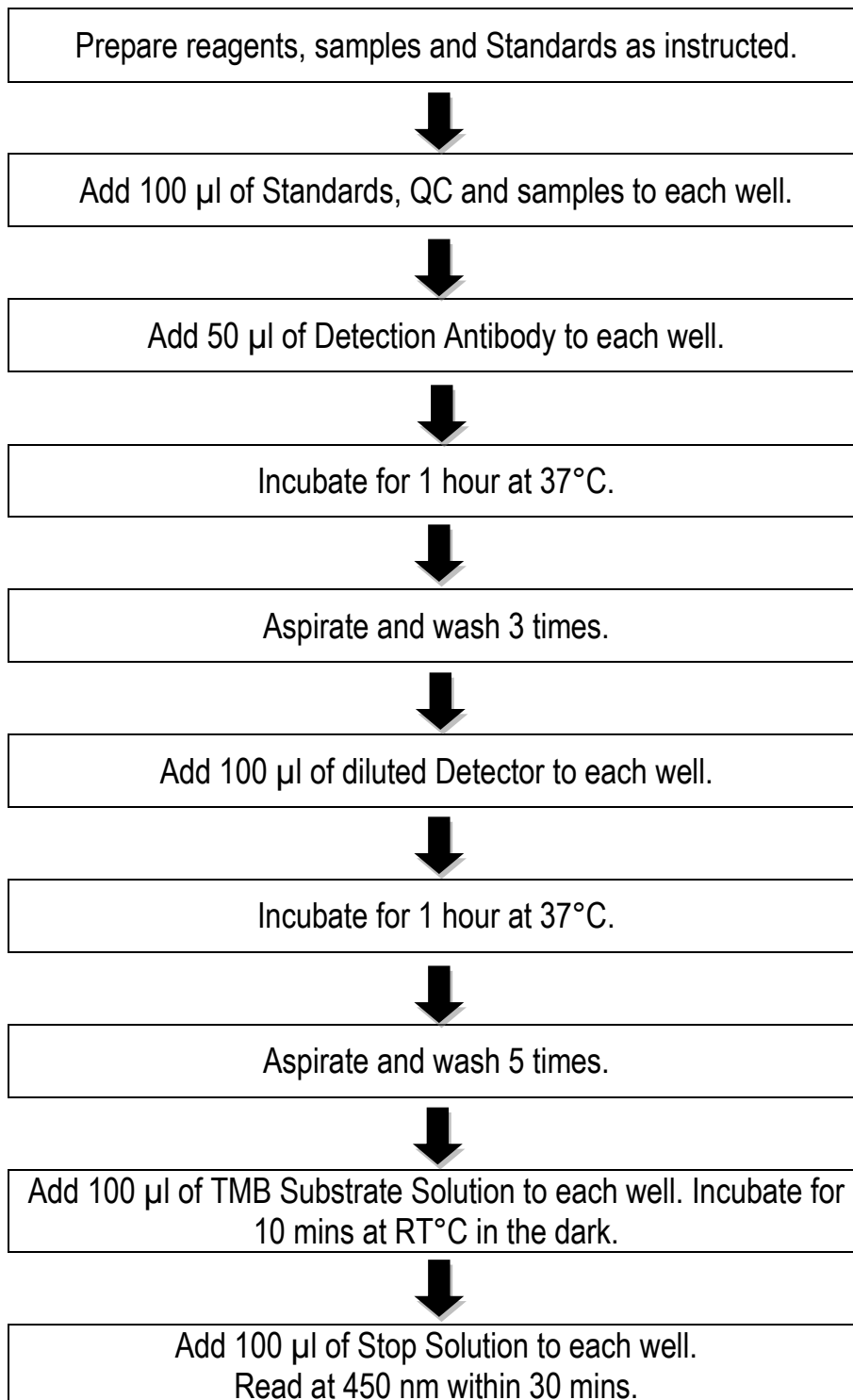


NOTES

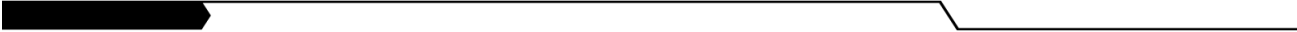

# BioVendor

Research  
and Diagnostic Products

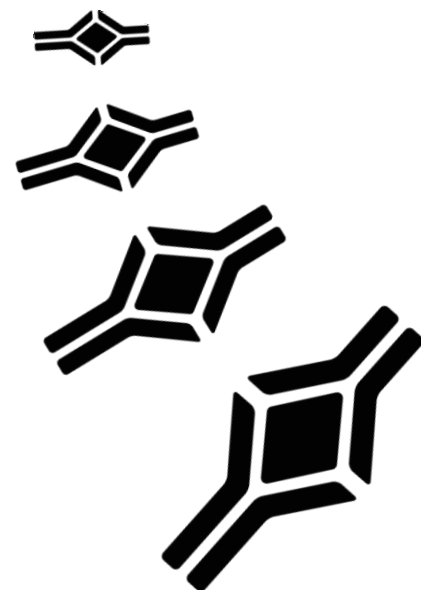

## HEADQUARTERS:

BioVendor - Laboratorní  
medicína a.s.

Karasek 1767/1

621 00 Brno  
CZECH REPUBLIC

Phone: +420-549-124-185  
Fax: +420-549-211-460

E-mail: [info@biovendor.com](mailto:info@biovendor.com)  
Web: [www.biovendor.com](http://www.biovendor.com)

EUROPEAN UNION:  
BioVendor GmbH

Im Neuenheimer Feld 583

D-69120 Heidelberg  
GERMANY

Phone: +49-6221-433-9100  
Fax: +49-6221-433-9111

E-mail: [infoEU@biovendor.com](mailto:infoEU@biovendor.com)

USA, CANADA AND MEXICO:  
BioVendor LLC

1463 Sand Hill Road  
Suite 227

Candler, NC 28715  
USA

Phone: +1-828-670-7807  
+1-800-404-7807  
Fax: +1-828-670-7809

E-mail: [infoUSA@biovendor.com](mailto:infoUSA@biovendor.com)

CHINA - Hong Kong Office:  
BioVendor Laboratories Ltd

Room 4008  
Hong Kong Plaza, No.188

Connaught Road West  
Hong Kong, CHINA

Phone: +852-2803-0523  
Fax: +852-2803-0525

E-mail: [infoHK@biovendor.com](mailto:infoHK@biovendor.com)

CHINA – Mainland Office:  
BioVendor Laboratories Ltd

Room 2917, 29/F  
R & F Ying Feng Plaza, No.2  
Huaqiang road

Pearl River New Town  
Guang Zhou, CHINA

Phone: +86-20-38065519  
Fax: +86-20-38065529

E-mail: [infoCN@biovendor.com](mailto:infoCN@biovendor.com)